



DERMATOPHYTE PROFESSIONAL TEST STRIP

TEST YOUR CLIENTS FOR TINEA UNGUIUM IN YOUR OWN CLINIC



SAVES EXPENSIVE COURSES



ONLY 10 MINUTES TIME



SAVES MONEY ON EACH TEST



COMPATIBLE WITH ANTIFUNGAL DRUGS

ARE YOU INTERESTED? DO NOT HESITATE TO CONTACT US.

WWW.HFLLABORATORIES.COM | WWW.FUNGICHECK.COM

HFL LABORATORIES BV | DE HOOGT 59 | 5175 AX LOON OP ZAND | NL T. +31 (0)13 4551749 | INFO@HFLLABORATORIES.COM

THE REVOLUTIONAIRY INNOVATION

Dermatophyte Test Strip FungiCheck Tinea Unguium



DERMATOPHYTE PROFESSIONAL TEST STRIP

TEST YOUR CLIENTS FOR TINEA UNGUIUM IN YOUR OWN CLINIC

HFL LABORATORIES | ESTABLISHED IN THE NETHERLANDS SINCE 2005

THE REVOLUTIONARY INNOVATION

FungiCheck is an immunochromatographic kit for professionals to detect Dermatophytes in a nail sample.

Dermatophytes are fungi that require keratin for growth. These fungi can cause infections of the skin, hair and the nails. Because Dermatophytes need keratin to grow, they are limited to hair, nails and the superficial skin. Dermatophytosis are referred to as "Tinea" infections.

The FungiCheck test strip is limited to the testing of the presence of Tinea Unguium (Onychomycosis) in the nails. (With the exception of Candida & yeasts)


Examples of Tinea unguium expressions:

- Lateral onychomycosis. A white or yellow strip on one side of the nail.
- Hyperkeratosis. Peelable layers beneath the nail.
- O Distal onycholysis. The end of the nail comes up. The free edge often crumbles.
- Superficial white onychomycosis. Flaky white patches and pits on the top of the nail plate.
- Proximal onychomycosis. Yellow spots appear in the Lunula.
- Complete destruction of the nail.

HOW DOES IT WORK?

The FungiCheck Test Kit includes 10 complete tests. Including 10 test strips, extraction liquid, 10 test tubes and 10 stirrers.

The test must always be performed by a professional with some experience in sterile testing. The test is very easy to use. With just 1mg nail sample and 10 minutes of your time you can already read the result of the test strip

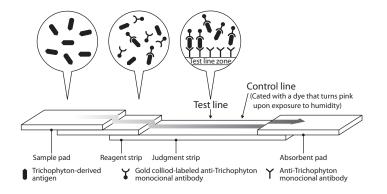

www.FungiCheck.com | www.HFLLaboratories.com

The test strip in the FungiCheck kit consists of a sample layer, a reagent strip, a decision strip and an absorbent pad. The extra control line assures you that your sample is appropriate absorbed by the test strip.



www.FungiCheck.com | www.HFLLaboratories.com